


Level Up Quality. Redefine Sharp.

Reusable & Disposable



NEW Slit Knife RIB&R

FEATHER® MICRO FEATHER Ophthalmic Scalpel

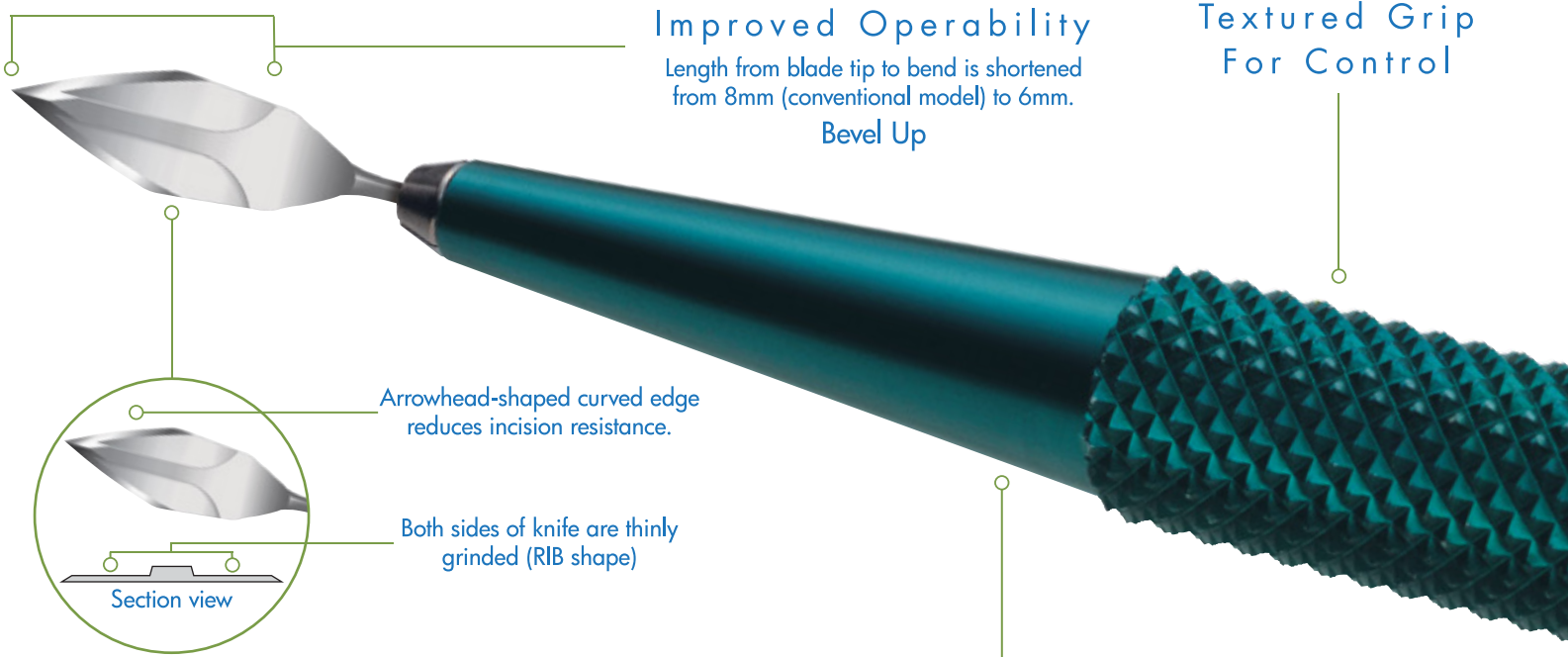


Improved Operability
Length from blade tip to bend is shortened from 8mm (conventional model) to 6mm.
Bevel Up

Textured Grip For Control

Arrowhead-shaped curved edge reduces incision resistance.

Both sides of knife are thinly grinded (RIB shape)

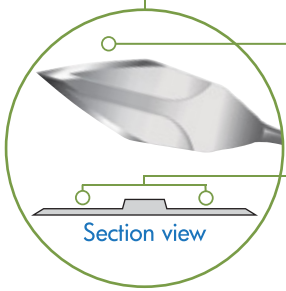

Lightweight, Aluminum Handle

Blade Width	2.2mm	2.4mm	2.75mm	2.8mm	3.0mm
Aluminum Handle	0622RC	0624RC	06275RC	0628RC	0630RC

MICRO FEATHER OPHTHALMIC SCALPEL Slit Knife <RIB & R> reduces incision resistance and allows smoother incision than ever before by its exclusively designed RIB-shaped knife and arrowhead-shaped curved edge.

Feather® is dedicated to manufacturing precision medical products that meet strict international quality standards.

Call to Order
(844) 820-8940

Learn how OASIS® can collaborate with your surgical team for **True Surgical Solutions.**

